SIEMENS

Behind-the-Ear Hearing Instruments

User Guide

usa.bestsoundtechnology.com

Life sounds brilliant.

Enhance your Siemens experience.

Learn about your hearing aids and accessories you may find helpful, access the latest information about hearing health, and receive promotional offers by joining the Siemens hearing health community at* usa.bestsoundtechnology.com/join or:

Learn more
about Siemens
products by
scanning the
QR code

Thank you for choosing our hearing instruments to accompany you through everyday life. Like anything new, it may take you a little time to become familiar with them.

This guide, along with support from your Hearing Care Professional, will help you understand the advantages and greater quality of life your hearing instruments offer.

If this Hearing Instrument is being purchased for an infant or young child, please confirm with your Hearing Care Professional that it is recommended for this purpose and to ensure appropriate supervision and use of the product.

Hearing Aid Information

Serial Number:

Left: _____ Right: _____

Battery Size:

Left: _____ Right: _____

Purchase Date:

Left: _____ Right: _____

*We respect your privacy. We will not sell, rent, or give your name or e-mail address to anyone.
At any time, you can select a link at the bottom of every email to unsubscribe or update profile.

Hearing Care Professional Info stamped here

Table of contents

| | |
|---|----|
| Getting started | |
| Your hearing instruments | 4 |
| Safety information | 4 |
| Feature card | 6 |
| Identification information | 6 |
| Identifying the left and right hearing instrument(s) | 7 |
| Changing the battery | 8 |
| Battery tips | 9 |
| Using the eCharger™ | 10 |
| Turning the hearing instrument on or off | 11 |
| Attaching the earmold | 11 |
| Inserting and removing hearing instrument with an earmold | 12 |
| Inserting and removing the hearing instrument with life tube or s-life tube | 13 |
| Operating your hearing instrument – Options and features | |
| Programmable and exchangeable controls | 15 |
| easyTek™ and miniTek® Wireless Enhancement, and easyPocket™ Remote Control | 17 |
| Adjusting the volume | 18 |
| Further adjustments (optional) | 18 |
| Locking the controls (optional) | 18 |
| Directional microphone system | 19 |
| Telecoil | 19 |

| | |
|---|----|
| AutoPhone®/AutoPhone magnet accessory | 20 |
| Audio input | 21 |
| Battery compartment door lock | 23 |
| Operating your hearing instrument – Programmable features | |
| Wireless functionality | 24 |
| Wireless information | 26 |
| Push button on/off | 28 |
| Power on delay | 28 |
| Low battery indicator | 28 |
| Volume control alert tones | 29 |
| Program change alert tones | 29 |
| Tinnitus therapy feature | 29 |
| Listening tips | 33 |
| Maintenance and care | 34 |
| Performance checklist | 36 |
| Health considerations | 37 |
| Important notices | 38 |
| Warranty and service | 40 |

Getting Started

Your hearing instruments

By wearing your hearing instruments you may once again enjoy the sounds of the environment and conversation with your family and friends. Your decision to purchase our hearing instruments reflects the importance you place on quality and reliability.

Your Behind-the-Ear (BTE) hearing instruments are designed for comfort, performance and durability. They should fit comfortably behind your ears and are attached to custom-made earmolds or flexible LifeTubes and LifeTips.

Please read this manual carefully to benefit from all the features of your new hearing system.

Safety information

WARNING

Choking hazard posed by small parts.

- Keep hearing instruments, batteries and accessories out of children's reach.
- If swallowed consult a physician or hospital immediately.
- If infants, small children or mentally disabled persons need to wear hearing instruments ensure adequate supervision.
- Ensure appropriate supervision when your child is wearing the hearing instruments.
- Check the integrity of the hearing instruments regularly.
- Ensure that your child does not detach the hearing instrument from the earmold.

- Consult your Hearing Care Professional if the housing is deformed.
- Keep the battery compartment locked when possible. Verify the proper function of the locking mechanism.

WARNING

Risk in explosive environments, e.g., mining areas

- Don't use your hearing instruments (and remote control) in areas where there is danger of explosions.

WARNING

- Risk of impairing the residual hearing of the user. Use only hearing instruments that have been fitted especially for your needs.
- Risk of electrical shock caused by damaged devices. Do not use obviously damaged devices and return them to your hearing health provider.
- Devices may interfere with electronic equipment. Verify your devices are turned off in areas where the use of electronics or wireless devices is restricted.

CAUTION

Risk of injury caused by incorrect use of Life Tip

- Always wear the Life Tube with the Life Tip.
- Make sure that the Life Tip is completely attached to the Life Tip adapter.
- Change the Life Tube or Life Tip when they are cracked.
- Carefully insert the Life Tube with the Life Tip into the ear canal. Be sure not to insert too deeply.

- Use the Life Tube and the Life Tip in the way described in this user guide.

If the Life Tip remains in the ear, have it removed by a medical professional.

► Feature card

Please refer to the feature card included with this manual. The feature card is located at the back of this user manual in a pocket.

The feature card contains important information about the features in your new hearing instruments. This includes battery size, and features associated with the operation of your hearing instruments.

More information about these features is located in the Operation section of this manual.

Identification information

Your hearing instruments have a serial number engraved on them. The location of the serial number will vary according to the

style of hearing instrument you have chosen. Record the serial number in the front of this manual for future reference (Fig. 1).

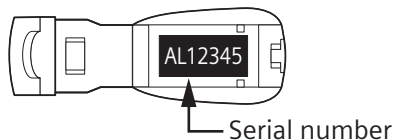
Fig. 1: Identification information

| Code (2nd Digit) | Year |
|---------------------|------|
| B, C | 2012 |
| D, E, F | 2013 |
| G, H | 2014 |
| L, M | 2015 |
| N, P | 2016 |

The year of manufacture is incorporated into the serial number. Please note the warranty is based upon the date of purchase, not the date of manufacture.

The serial number may look similar to this sample: AL12345. The second letter of this serial number refers to the year of manufacture of the hearing instrument. The letters L and M indicate 2015 year of manufacture, while N and P indicate 2016. Your Hearing Care Professional can answer any questions you may have about the identifying code on your hearing instrument.

Identifying the left and right hearing instrument(s)

The left ear hearing instrument may have a blue typeplate or a blue indicator inside the battery compartment. The right ear hearing instrument may have a red typeplate or a red indicator inside the battery compartment (Fig. 2).

In addition, on BTEs, the left ear hearing instrument has a blue number on the LifeTip Adapter, and the right ear hearing instrument has a red number on the LifeTip Adapter (Fig. 3).

Fig. 2: Indicator

Fig. 3: Location of right and left designator on LifeTip adapter for Siemens Life™ Product

Changing the battery

Refer to the feature card or check with your Hearing Care Professional to determine the battery size for your hearing instrument.

To insert the battery into the hearing instrument, use your fingernail to swivel open the door of the battery compartment. Place the proper size battery in the battery compartment so that the “+” symbol on the battery coincides with the “+” mark on the battery compartment door.

Gently close the battery compartment. Do not force the door shut. If it does not close easily, check to see if the battery was inserted properly. When the battery door is completely closed, your hearing instrument is ready for operation (Fig. 4, 4.1, and 4.2).

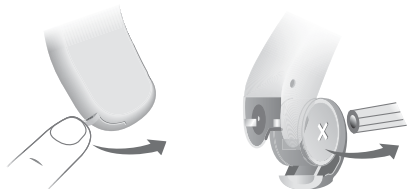
Fig. 4: Inserting the battery

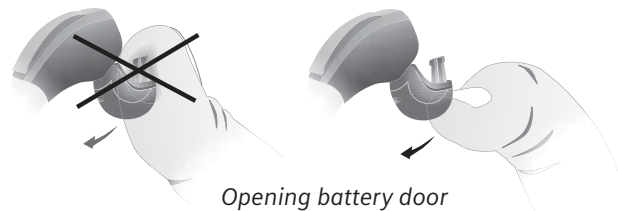
Fig. 4.1: Inserting the battery


Your hearing instrument may be programmed to inform you when the battery should be replaced.

Refer to the low battery indicator information in the Operation section (page 23) of this manual for further information.


Fig. 4.2: Changing the battery on Aquaris

Opening battery door

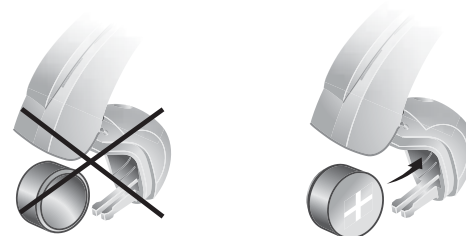Inserting battery

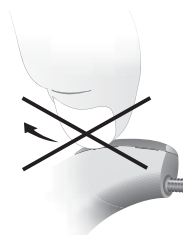The microphone cover on Aquaris is not detachable

Battery tips

- If a battery is accidentally swallowed, seek medical attention immediately, or call the National Battery Hotline collect at (202) 625-3333.

- Always use the recommended battery size for your hearing instrument. The battery size information for your hearing instrument is located on the feature card.
- Remove the battery from your instruments when not in use. This prevents possible damage from expansion or leakage of the battery, if the instruments are not fully turned off.
- Always carry spare batteries in the event that your current battery goes dead. Place each spare battery in an individual, non-metallic container to prevent contact with other batteries or metal objects, such as keys or coins.
- Follow your local recycling rules for battery disposal.

Using the eCharger

If you have the optional charger, refer to the charger instruction manual for operating instructions. Use only the rechargeable batteries delivered with the charger or available through your Hearing Care Professional. The lifetime of a rechargeable battery is approximately one year, depending on individual use.

NOTE: Do not attempt to charge standard, non-rechargeable batteries!

NOTE: The hearing instrument must be turned on after removing from the charger (see page 10 Turning the Hearing Instrument on or off).

Fig. 4.3: Charger

Turning the hearing instrument on or off

To turn the hearing instrument on, insert a battery into the hearing instrument and completely close the door.

To turn the hearing instrument off, open the battery door to the first catch.

Your hearing instruments may also be programmed to allow you to turn them on and off using the program button. Refer to the Operation section (page 20) for further information.

You can also use an optional remote control to turn the hearing instruments on and off.

Some BTE's may also have a LED light (Fig.5) near the battery door. The LED (optional) shows whether the hearing instruments are on. LED can be covered or programmed on/off

- The LED shines continuously: the hearing instruments are switched on and sufficient battery power is available.
- The LED is off: The hearing instruments are switched off or the battery must be replaced.

Fig. 5: LED light

Attaching the earmold (where provided)

If your earmold detaches, gently insert the earhook of the hearing instrument into the tubing of the earmold supplied by your Hearing Care Professional (Fig. 6).

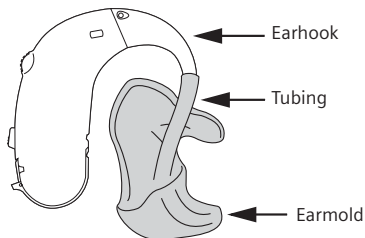
Fig. 6: Attaching the earmold

Inserting and removing hearing instruments with an earmold

Hold the earmold between your thumb and index finger. Gently work the earmold into its proper position by adjusting it until it is firmly seated in your ear. Then lift the hearing instrument over the top of your ear and adjust it to fit behind your ear (Fig. 7).

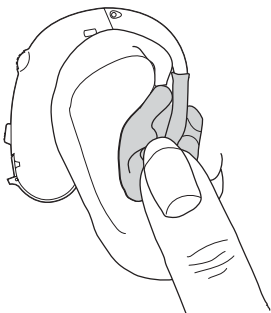
Fig. 7: Inserting the hearing instrument

Lightly push on the back of the ear to help release the earmold and gently pull it out (Fig. 8).

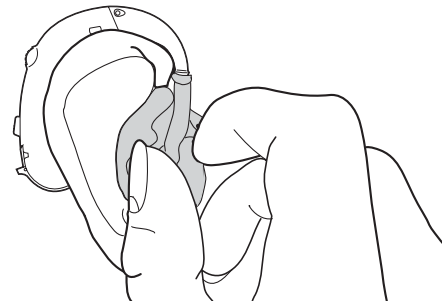
Fig. 8: Removing the hearing instrument

Inserting and removing hearing instruments with LifeTube or S-Life Tube

Hold the LifeTip between your thumb and index finger. Gently work the LifeTip into its proper position in the ear canal by adjusting it until it is firmly seated. Then lift the hearing instrument over the top of your ear and adjust it to fit behind your ear. Finally, place the retention strand into the bottom of the "bowl" of your ear (Fig. 9).

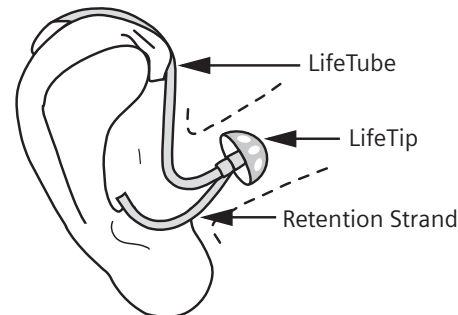
Fig. 9: Inserting the hearing instrument

To remove the instrument, gently pull the LifeTip out of your ear and carefully lift the hearing instrument off your ear. Do not pull or tug on the LifeTube to remove the Tip or instrument, or the LifeTube may become damaged (Fig. 10).


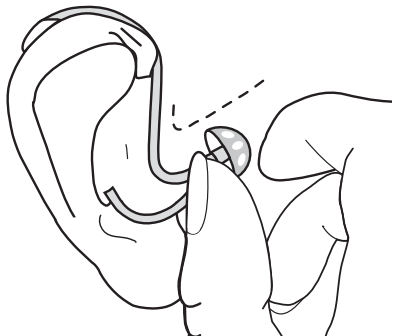
Fig. 10: Removing the hearing instrument

Operating Your Hearing Instrument – Options and Features

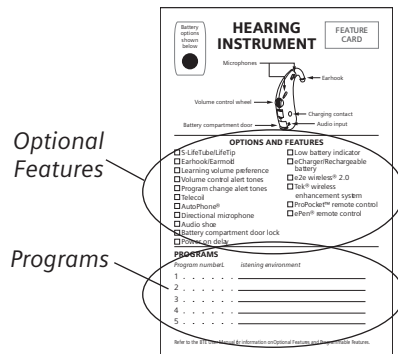Optional Features

Programs

Refer to the feature card included with this manual to determine the optional features your hearing instrument is equipped with (Fig. 11):

Fig. 11: FEATURE CARD

Options and Features section circled

Programmable and Exchangeable Controls

With the controls you can, for example, adjust the volume or switch hearing programs. Your hearing instruments either have a push button or a rocker switch – or no control at all.

Your Hearing Care Professional has assigned your desired functions to the controls, and has marked their functions on your feature card.

Rocker Switch – If your rocker switch is configured as a volume control pushing it up makes the volume louder. Pushing it down makes the volume softer.

Your hearing instrument may have several memories, also called hearing programs.

If your Rocker Switch is configured as a Program Button then

pushing it up switches to the next program. Pushing it down switches to the previous program.

If the Rocker Switch is configured as a Volume Control or Push Button:

Press 1 second to change Volume

Press 2 seconds to change Program

Press 3 seconds to turn on/off

Tinnitus Therapy Feature – You need to be in the Tinnitus Therapy Feature program. Press the rocker switch briefly to increase or decrease the tinnitus therapy feature level.

When initially turned on (after the battery compartment door is closed), the microphone is ON and ready for general use. The default is Program 1. Your Hearing Care Professional will provide you with the number of programs most appropriate for your needs. Each time you press the program button or rocker switch on the hearing instrument (if available) or optional remote control the instrument changes to the next program. The program change may be confirmed by alerting tones or beeps that correlate to the program you are in. That means if you select Program 3, then the hearing instrument beeps 3 times.

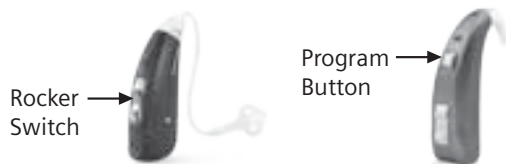
Fig. 12: Changing the program with the program button or rocker switch

Your Hearing Care Professional may list your personal programs on the feature card (Fig. 11).

easyTek™ and miniTek® Wireless Enhancement, and easyPocket™ Remote Control

You may have received the optional easyTek, miniTek, or easyPocket (Fig. 14) with your hearing instrument. miniTek or easyPocket can be used to change the program and volume of your hearing system. In addition, other adjustments can also be made to sound quality.

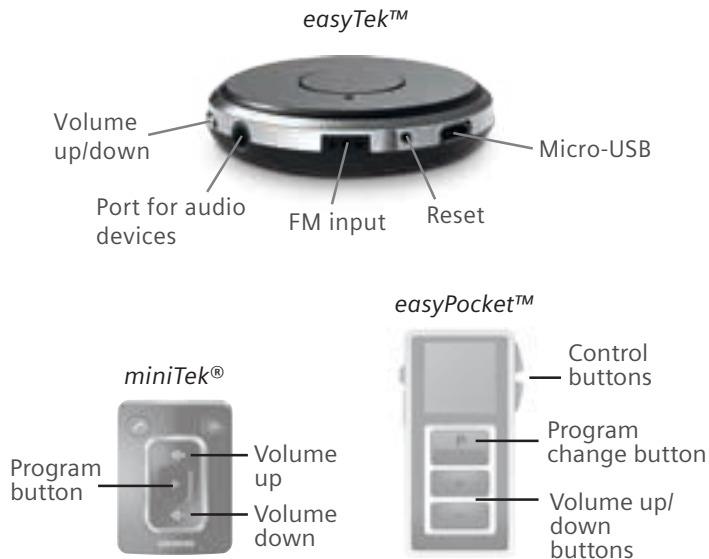
Fig. 14: easyTek, miniTek, and easyPocket remote controls

For detailed information and instructions refer to the easyTek User Guide, miniTek User Manual or the Remote Control User Manual for easyPocket or ask your Hearing Care Professional.

Adjusting the volume

Your hearing instruments automatically adjust the volume to the listening situation.

If you prefer manual volume adjustment, press the push button or the rocker switch, or use a remote control.

An optional signal can indicate the volume change.

Further adjustments (optional)

The controls of your hearing instruments can also be used to change the sound balance or the tinnitus therapy feature level. Sound balance allows you to adjust the bass or treble to your preferences.

Locking the controls (optional)

To prevent any accidental operation you can lock the controls. When the controls are locked, both push button and rocker switch are disabled.

- To lock the controls, press the top of the rocker switch and hold it. Simultaneously press the push button on the same hearing instrument for 3 seconds.
- To unlock the controls, press the bottom of the rocker switch and hold it. Simultaneously press the push button on the same hearing instrument for 3 seconds.

Whether this function is configured for your hearing instruments is listed in section "Controls".

Directional microphone system

According to your hearing needs, your hearing instrument may be adjusted to focus on sounds in front of you (directional mode) for when you are having a face-to-face conversation in a noisy environment or behind you when the speech you wish to hear is coming from the back such as in a car. The hearing system may also be adjusted for when you are in more quiet environments and want to hear the sounds around you (omnidirectional). Your Hearing Care Professional may program the operation of the directional microphone system to function automatically or in conjunction with the program button or remote control based on your particular hearing requirements.

Telecoil

Your hearing instrument may be equipped with a telecoil for telephone use or use in looped environments. It may be programmed to work in a specific listening program in conjunction with the program button. In the telecoil mode, place the telephone receiver close to your hearing instrument for best reception. By trying various positions of the telephone receiver, you will discover what works best for you. Your Hearing Care Professional can advise you on the proper use of the telecoil.

Note: In order for the telecoil to function properly, your telephone must be "hearing aid compatible". Check your telephone's instruction manual to find out whether your phone has this capability.

AutoPhone®

The AutoPhone option allows the telecoil to automatically engage when a "hearing aid compatible" telephone handset is brought close to the hearing instrument.

To operate, bring the telephone handset close to the hearing instrument. The telecoil will automatically engage to the presence of the electromagnetic signal from the telephone. When you are finished speaking on the phone, move the handset away from the hearing instrument and the microphone mode will re-establish itself.

AutoPhone magnet accessory

There is an AutoPhone magnet accessory available for your telephone that can be used with the AutoPhone feature of your hearing instrument. Ask your Hearing Care Professional about this accessory and if it is appropriate for your individual hearing care needs.

You must use the original magnet provided with your hearing instrument. Using other magnets may cause irreversible damage to the electronics of your hearing instrument.

Note: When magnets are properly cared for, they are safe to use and can last a lifetime. However, the following guidelines should be followed while using these magnets.

- Magnets can pose risks to electronics and magnetic storage media when placed in close proximity.
- Be careful not to place magnets near computers, color monitors, television sets, floppy disks, video/audio tape, and other electronic equipment.
- Magnets may interfere with the operation of life support systems, such as pacemakers, magnetic valves, etc.
- Make certain to maintain a distance of at least 4 inches between the magnet and the area of your support system.
- Keep magnets out of the reach of children or pets.
- If accidentally swallowed, seek medical attention.

Audio input

Your hearing instruments may be equipped with audio input capability. This enables you to connect your hearing instruments to other audio sources such as an FM system, tape recorder, television, or infrared system. It can be useful for adults in difficult listening situations and for children in classrooms. Because the sound is coming directly from the sound source to the hearing instrument, the signal is clearer and unaffected by poor room acoustics. Ask your Hearing Care Professional if this feature may be appropriate for your needs. There are two methods for attaching and removing the audio input shoe. Refer to the next page for the method you should follow.

WARNING

Risk of injuries caused by electric shock!

Connect audio input only to devices conform to IEC 60065 (IEC-Standard for audio, video and other electronic devices).

Attaching the audio shoe

1. Open the battery compartment.
2. Pull the cover out. The shape of the cover depends on your type of hearing instrument.

3. Slide the audio shoe and ensure that it clicks noticeably in place.

4. Close the battery compartment.

Removing the audio shoe

When the audio shoe is not attached, you can install a cover to prevent dirt and dust from entering the housing.

1. Open the battery compartment.

2. Drag the slider on the audio shoe back and at the same time move the audio shoe downwards to remove it from the device.

3. Slide the cover in the direction of the arrow into position. The shape of the cover depends on your type of hearing instrument.

4. Close the battery compartment.

Battery compartment door lock

Your hearing instruments can be modified to provide a child-safe battery compartment, which can be locked (child lock). The child lock has to be activated by the Hearing Care Professional.

With the child lock enabled, the battery compartment will automatically lock when closed. This prevents e.g. a child from removing and swallowing the battery.

To open the locked battery compartment, you need a tool. The tool is included with delivery of any child-lock battery compartment set.

Opening a locked battery compartment:

- Use the tool to push the white locking bar to the opposite side.
- Hold it in place and at the same time open the battery compartment.

Locking the battery compartment:

- Insert the tip of the tool into the small opening.
- Push / Lever the small plastic part out of the battery compartment.

Wireless functionality (if applicable)

Your hearing instrument may be equipped with wireless functionality. This provides synchronization between two hearing instruments. If you change the program or volume on one instrument, the other is automatically changed. In addition, wireless technology automatically adjusts the system to change both hearing instruments simultaneously to adapt to different listening situations. The wireless functionality also enables use of an optional remote control.

Note: The performance of the e2e® wireless 2.0 could be temporarily affected when in close proximity to electromagnetic noise sources, such as a computer monitor or television. Move away from the noise source if you experience difficulty.

Technical information for e2e 2.0 wireless function:

Operating frequency : $F_c = 3.28 \text{ MHz}$; Maximum field strength: $-7 \text{ dB } \mu\text{A/m}$ at 10m meters.

Changes or modifications not expressly approved by the party responsible for compliance could void the users authority to operate the equipment.

This Class B digital apparatus complies with Part 15 of FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause interference, and (2) this device must accept

any interference, including interference that may cause undesired operation of the device.

FCC: SGI-WL002BTE.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules.

These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

Wireless Information

The following table summarizes the technical details of the wireless technology:

| | Wireless Hearing Aids |
|--|--|
| Antenna Type | Inductive Antenna |
| Antenna Dimensions | Ø1.9 mm, L – 6.5 mm |
| Modulation | PSK (Phase Shift Key) |
| Magnetic Field Strength (at 10 m distance) | 0.09 A/m (1 cm ² coil; average) |
| Output Power (EIRP**) | 42 µW |
| Range | <1 m remote to hearing aid <30 cm between aids or aid to remote |
| Center Frequency | 3.28 MHz |
| Channel | Single Channel Radio |
| Bandwidth | 67 kHz |
| Data-Rate | 655 kbit/sec (raw channel capacity) |
| Data Flow | Simplex or semi-duplex capability |
| Protocol | Random access – no collision avoidance |
| S.A.R.* | 5 nW/kg (BTE), 13 nW/kg (ITE) |

*Specific Absorption Rating (S.A.R.) for continuous transmission (100% duty cycle) based on 10g ICNIRP testing.

EMI/EMC COMPLIANCE

Wireless hearing aids comply with the following EMC/EMI standards:

| Standard | Test Type | Note |
|---------------------------|--|---|
| 47 CFR Part 15, Subpart C | RF emissions | U.S. FCC requirements for intentional radiators. |
| EN 300 330-1/2 | RF emissions including spurious emission | EMC and radio spectrum matters for short range devices in the frequency range 9 kHz – 25 MHz. |
| IEC 60601-1-2 | EMC emission immunity, RF and ESD | Medical electrical equipment. General requirements for basic safety and essential performance. Electromagnetic compatibility. |
| EN 301 489-1/3 | Immunity, RF and ESD | Standard for low power transmitters in the frequency range 9 kHz – 40 GHz. |
| IEC 60118-13 | Immunity RF near field immunity test | International product standard for hearing aids to ensure adequate immunity to radio interference from mobile telephones. |
| ANSI C63.19 | Immunity RF near field immunity test | American National Standard method of measurement of compatibility between wireless communication devices and hearing aids. |
| ANSI/AAMI PC69:2007 | EMC immunity | Implantable medical device EMC immunity. |
| EN 45502-2-1 | EMC immunity | Particular requirements for pacemakers. |

WIRELESS SECURITY MEASURES

Wireless signal security is assured through device system design that includes:

- A built-in pairing table which specifies valid and legitimate pairing among units.
- A proprietary communication protocol which checks the package numbers during each transmission.
- A Cyclic Redundancy Check (CRC) to check data validity and correct errors.

Push button on/off

Your hearing instruments may be programmed to allow you to turn them on and off using the program button. To turn the instrument off using the program button, press and hold the button for 3–5 seconds. To turn it back on, press and hold the program button for 3–5 seconds. Please note that after having turned the hearing instrument off using the program button, you must wait approximately 10 seconds for the system to reset before turning it back on again by using the program button.

Power on delay

Your hearing instruments may be equipped with an advanced power-on delay feature. This feature delays sound from entering your hearing instruments when they are first turned on so that you may insert the instruments into your ears without experiencing unpleasant feedback (whistling).

Low battery indicator

Your hearing instrument may be programmed to inform you when the battery should be replaced for best performance. You may hear an intermittent series of beeps and/or experience low

amplification power or volume. Change the battery immediately if either of these should occur.

Volume control alert tones

Your hearing instrument may be programmed to inform you of changes to volume. You may hear soft beeps each time you change the volume and a three beep warning when the volume limit has been reached. Once the maximum or minimum volume has been reached, further changes in that direction will have no effect.

Program change alert tones

Your hearing instrument may be programmed to inform you of changes to the program. The program change may be confirmed by alert tones or beeps that correlate to the program you are in. That means if you select Program 3, the hearing instrument will beep 3 times.

Tinnitus Therapy Feature

TINNITUS THERAPY FEATURE:

Your hearing instrument may be equipped with a tinnitus therapy feature. Many tinnitus patients also suffer from some degree of hearing loss. Many hearing instruments can be used alone or in combination with the tinnitus therapy feature. The tinnitus therapy feature is fixed to a broadband noise that can be adjusted by your Hearing Care Professional for your tinnitus therapy. The feature may provide temporary relief of your tinnitus.

Prescription Use Only

Federal law restricts this device to sale by or on the order of a doctor or hearing care professional licensed to dispense hearing instruments in your state. The use of any sound generating tinnitus therapy device should be only on the advice and in consultation with your hearing care professional. Your hearing care professional will properly diagnose and fit the device to your personal needs and requirements. This should include its use in a prescribed tinnitus treatment program. Your hearing care professional will also be able to offer the appropriate follow-up care. It is important that you follow your hearing care professional's advice and direction regarding such care.

WARNING

Risk of further impairment to the user's hearing health.

There are some potential concerns associated with the use of any sound generating tinnitus therapy device. Discontinue use and seek medical evaluation if any of the following conditions occur:

- Chronic skin irritation on, near, or around the site of device placement.
- Unusual side effects (e.g., dizziness, nausea, headaches, heart palpitations).
- Perceived decrease in auditory function (e.g., decreased loudness, speech not as clear).

WARNING

Risk of further impairment to the user's hearing health.

The volume of the tinnitus therapy feature can be set to a level which could lead to permanent hearing damage when used for a prolonged period of time.

The tinnitus therapy feature should never be used at uncomfortable levels.

FOR HEARING CARE PROFESSIONALS

Indications for Use

The tinnitus therapy feature is a tool to generate sounds to be used in a tinnitus management program to relieve patients suffering from tinnitus. The target population is primarily the adult population over 21 years of age. The tinnitus therapy feature is targeted for healthcare professionals treating patients suffering from tinnitus, as well as conventional hearing disorders. The fitting of the tinnitus therapy feature should be done by a hearing care professional participating in a tinnitus management program.

Device Description

The tinnitus therapy feature is a software function that generates sound which is programmed into a hearing instrument. Depending on the type of hearing instrument, the hearing instrument may be used in up to three modes of operation: as a hearing instrument, as a tinnitus treatment device, or as a hearing instrument and tinnitus treatment device.

When enabled, the tinnitus therapy feature generates the sound and allows a patient's hearing care professional to design and program appropriate settings for an individually prescribed sound treatment plan. The treatment plan should be used in a tinnitus management program for relief of tinnitus.

The tinnitus therapy feature generates a broadband noise signal that varies in frequency and amplitude. These characteristics are adjustable by the hearing care professional and are specific to the prescribed therapy designed by the professional for the patient's needs and comfort. The patient may have some control of the level or volume of the signal and the patient should discuss this adjustment as well as his or her comfort level and sound of the signal with their hearing care professional.

WARNING

Risk of further impairment to the user's hearing health.

The volume of the tinnitus therapy signal can be set to a level which could lead to permanent hearing damage when used for a prolonged period of time.

Should the tinnitus therapy signal be set to such a level in the hearing instrument, advise the user of the maximum amount of time per day he or she should use the tinnitus therapy signal. For example, occupational safety guidelines restrict continuous noise exposure of 80dBA SPL to 8 hours per day. The tinnitus therapy signal should never be used at uncomfortable levels.

WARNING

Risk of further impairment to the user's hearing health.

A hearing care professional should advise a prospective tinnitus therapy feature user to consult promptly with a licensed physician (preferably an ear specialist) before using the tinnitus therapy feature if the hearing care professional determines through inquiry, actual observation, or review or any other available information concerning the prospective user that the prospective user has any of the following conditions:

- i. Visible congenital or traumatic deformity of the ear.
- ii. History of active drainage from the ear within the previous 90 days.
- iii. History of sudden or rapidly progressive hearing loss within the previous 90 days.
- iv. Acute or chronic dizziness.
- v. Unilateral hearing loss of sudden or recent onset within the previous 90 days.

Listening Tips

- Give yourself time to become accustomed to your new hearing instruments. During this period, you may prefer to wear your hearing instruments for only a portion of the day, and then gradually increase your usage throughout your normal routine each day.
- You may find that a quiet environment is more comfortable in the beginning of the adjustment period. Gradually introduce new listening environments as your ears become accustomed to hearing instrument amplification.
- If you are a new hearing instrument user, be aware that there are sounds, both pleasant and unpleasant, that you may not have heard in a long time.
- Ask the people with whom you are conversing to attract your attention before speaking to you. If you don't understand what someone has said, ask the speaker to slow down or say it in another way, but not to shout.
- In noisy or difficult environments, you may find it easier to understand speech by facing the speaker. Focus on reading lips and observing facial expressions and body language. Be aware that in some situations, it may be too noisy to understand what is being said.
- Encourage one person to speak at a time when you are engaging in group conversation.
- Listen to radio and television programs as you normally would.
- Experiment finding the best place to sit for optimum hearing in public places like the theater, school auditoriums, religious services, etc. If possible, choose a seat that is far away from any source of distracting noise, such as a fan or an open window, and where you have a good view of the person speaking.

Maintenance and Care

Hearing instruments should be cleaned daily. You can wipe them with a tissue or a soft cloth. Do not use water, alcohol, benzene, or any other liquid or household cleaning agent to clean your instruments as it may cause damage to them.

Always store your hearing instruments (and earmolds, LifeTube, or LifeTip) inside the case. A soft carrying pouch is also included for your convenience when traveling.

Remove the batteries to prevent damage from corrosion when storing hearing instruments for an extended period of time.

Store the instruments in a dry area, preferably at room temperature. If you perspire heavily or live in an area of high humidity, you may find it beneficial to use a DRI-AID™ kit when the instruments are not in use. This product can be purchased from your Hearing Care Professional.

Note: The charger can be used instead of the DRI-AID kit. Refer to the Maintenance and Care section of the charger instruction manual for further information.

Note: Your hearing instrument is sensitive to extreme heat, high humidity, strong magnetic fields (>0.1T), X-rays and mechanical stress.

- **Never immerse your hearing instruments in water. This does not apply to micon hearing instruments.** micon instruments are IP67 rated per IEC 60529 standard and are water-resistant devices that can be completely submerged in water up to 3 feet for 30 minutes. Dust will not interfere with the satisfactory operation of the device. Aquaris™ exception – see below for Aquaris instruments recommendations*.
- Do not expose your hearing instrument to extreme temperature

or high humidity. Aquaris™ exception. See below for Aquaris instruments recommendations*

- Do not leave it in direct sunlight.
- Do not wear it in the shower or when you apply makeup, perfume, aftershave, hairspray or suntan lotion. Aquaris exception. See below for Aquaris instruments recommendations*
- Do not wear your hearing instruments when you are exposed to short-waves, a strong magnetic field, a high frequency field or X-rays.
- Do not place your hearing instruments in a microwave oven.

*Aquaris™ Instruments Recommendations

Aquaris instruments achieved International Protection Rating (IP68) per the International Electrotechnical Commission (IEC) 60529 Standard. Your instrument can be completely submerged in water up to 3 feet with no damage to the instrument, and dust will not interfere with the operation of the device.

Special recommendations for water activities:

- When you are swimming, use an earhook with a custom-made, unvented, soft ear mold.
- Remove hearing instruments before diving and playing water sports that cause your hearing instruments to be constantly immersed in water.
- Water will not harm Aquaris, however, you should avoid getting water in your ear canal as the only entry point for water into Aquaris is through the receiver opening. Avoid applying shampoo and suntan or other lotions directly on Aquaris as it can impact its performance. Afterwards, dry Aquaris by gently blowing remaining water from under the microphone cover.

*Achieved IP68 rating per IEC 60529 standard. The device can be completely submerged with no water seeping inside, and no damage to the instrument during continuous immersion in water up to 3 feet. Also dust will not interfere with the satisfactory operation of the device Requires appropriate earmold for submersion.

Performance Checklist

| Problem | Cause | Possible Solution |
|---|---|--|
| Sound is weak | Weak or dead battery | Insert new battery |
| | Volume too low | Clean instrument, earmold, or change LifeTip |
| | Earmold or LifeTip clogged with wax or debris | Clean instrument, earmold or change LifeTip |
| | Microphone clogged | Wipe hearing instrument with cloth |
| Instrument emits whistling sound | Improper seating in ear | Try reinserting earmold or LifeTip until it fits securely |
| | Volume too high | Lower the volume |
| | Clogged with wax or excessive wax in ears | Clean earmold or replace LifeTip OR contact your Hearing Care Professional |
| Sound is weak | Weak or dead battery | Insert new battery |
| | Volume too low | Clean instrument, earmold, or change LifeTip |
| | Earmold or LifeTip clogged with wax or debris | Clean instrument, earmold or change LifeTip |
| | Microphone clogged | Wipe hearing instrument with cloth |
| Instrument emits whistling sound | Improper seating in ear | Try reinserting earmold or LifeTip until it fits securely |
| | Volume too high | Lower the volume |
| | Clogged with wax or excessive wax in ears | Clean earmold or replace LifeTip OR contact your Hearing Care Professional |

| | | |
|--|---|---|
| Hearing instrument emits signal tones | Battery is drained | Replace drained battery |
| Sound is distorted | Volume too high | Lower the volume |
| | Weak battery | Replace battery |
| Instrument does not work | Instrument not turned on | Turn instrument on |
| | Battery compartment not completely closed | Gently close the battery compartment completely |
| | Dead battery | Replace battery |
| | Battery polarity reversed | Replace battery |
| | Power-on delay is active | Wait 15 seconds and check again |

If the solutions in this checklist do not address a particular issue you are having with your hearing instruments, consult your Hearing Care Professional.

Health Considerations

If soreness or skin irritation develops, discontinue wearing your hearing instrument, and bring the instrument and earmold to your Hearing Care Professional. Minor fit adjustments or earmold modification can often correct this condition. If soreness persists, discontinue wearing the hearing instrument and see your physician. If excessive earwax accumulates when wearing your hearing instruments, consult your Hearing Care Professional.

Important Notices

Children with hearing loss

In addition to seeing a physician for a medical evaluation, a child with a hearing loss should be directed to an Audiologist, since hearing loss may cause problems in language development and the educational and social growth of a child.

An Audiologist is qualified by training and experience to assist in the evaluation and rehabilitation of a child with a hearing loss.

Important notice for prospective hearing instrument wearers

Good health practice requires that a person with a hearing loss have a medical evaluation by a licensed physician (preferably a physician who specializes in diseases of the ear) before purchasing hearing instruments. Licensed physicians who specialize in diseases of the ear are often referred to as Otolaryngologists, Otologists, or Otorhinolaryngologists.

The purpose of a medical evaluation is to assure that all medically treatable conditions that may affect hearing are identified and treated before hearing instruments are purchased.

Following the medical evaluation, the physician will give you a written statement that your hearing loss has been medically evaluated and that you may be considered a candidate for hearing instruments.

The physician will refer you to a Hearing Care Professional, as appropriate, for a hearing instrument evaluation.

The Hearing Care Professional will conduct a hearing instrument evaluation to assess your ability to hear with and without hearing instruments. The hearing instrument evaluation will enable the

Hearing Care Professional to select and fit hearing instruments for your individual needs.

If you have reservations about your ability to adapt to amplification, you should inquire about the availability of a trial-rental or purchase-option program. Many Hearing Care Professionals now offer programs that permit you to wear a hearing instrument for a period of time for a nominal fee after which you may decide if you want to purchase the hearing instrument.

Only those fully informed adults who have either obtained medical evaluation from a licensed physician or signed a waiver statement declining the medical evaluation for religious or personal beliefs that preclude consultation with a physician may purchase hearing instruments. The exercise of such a waiver is not in your best health interest and its use is strongly discouraged.

Warning to Hearing Care Professionals

An Hearing Care Professional should advise a prospective hearing instrument user to consult promptly with a licensed physician (preferably an ear specialist) before dispensing a hearing instrument if the Hearing Care Professional determines through inquiry, actual observation, or review of any other available information concerning the prospective user that the prospective user has any of the following conditions:

- Visible congenital or traumatic deformity of the ear.
- Acute or chronic dizziness.
- History of active drainage from the ear within the previous 90 days.
- History of sudden or rapidly progressive hearing loss within the previous 90 days.
- Unilateral hearing loss of sudden or recent onset within the previous 90 days.

- Audiometric air-bone gap equal to or greater than 15 decibels (dB) at 500 Hertz (Hz), 1000 Hz, and 2000 Hz.
- Visible evidence of significant cerumen accumulation or a foreign body in the ear canal.
- Pain or discomfort in the ear.

Special care should be exercised in selecting and fitting a hearing instrument whose maximum sound pressure level exceeds 132 decibels (dB), because there may be risk of impairing the remaining hearing of the hearing instrument user.

A hearing instrument will not restore normal hearing and will not prevent or improve a hearing impairment resulting from organic conditions.

The use of hearing instruments is only part of hearing rehabilitation and may need to be supplemented by auditory training and instructions in lip reading. Consistent use of amplification is recommended. In most cases, infrequent use of hearing instruments does not permit the wearer to attain the full benefit of amplification.

Warranty and Service

Your BTE hearing instrument, with the exception of the battery, is covered by a comprehensive warranty. All covered instrument parts received for warranty service at an authorized Service Center will be repaired or replaced with new or reconditioned components, without charge, to meet the performance specifications for that model.

This warranty does not cover malfunctions due to unusual wear and tear or mistreatment of the instrument such as physical

shock, excessive wax build-up, or tampering with the instrument, any of which voids all warranties. Your Hearing Care Professional may charge a service fee for processing warranty service.

Warranty service must only be performed by an authorized Service Center. Service performed by unauthorized service depots voids this warranty and repairs so necessitated will be done on a parts and labor cost basis.

Please refer to the warranty card included with your hearing instruments for warranty period effective dates.

Your BTE hearing instruments may have additional loss and damage coverage. Please consult your Hearing Care Professional to determine if this is applicable to your hearing instruments.

Procedure

In case of malfunction or loss, take your hearing instrument and/or proof of purchase to the Hearing Care Professional from whom the unit was purchased. If factory service or replacement is needed, your Hearing Care Professional will forward the unit to an authorized Service Center.

Extended warranty programs

We offer extended warranty and service programs for your hearing instruments. Please consult with your Hearing Care Professional for more details on these programs.

IMPORTANT

This warranty gives you specific legal rights and you may have other rights which may vary from location to location. For more information, contact your local Department of Consumer Affairs.

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usa.bestsoundtechnology.com/join o:

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Information in this brochure is subject to change without notice. Copyright © 2015 Siemens AG. All rights reserved.
La información que contiene este folleto está sujeta a cambios sin aviso. Copyright © 2015 Siemens AG. Todos los derechos reservados.
under Trademark License of Siemens AG
Sivantos, Inc.